


Hyperbaric Chamber Specifications



Four double layered polycarbonate windows allow for ample light to enter the chamber. The two face windows are offset so an attendant can observe the chamber patient from either side without having to lean over the chamber.

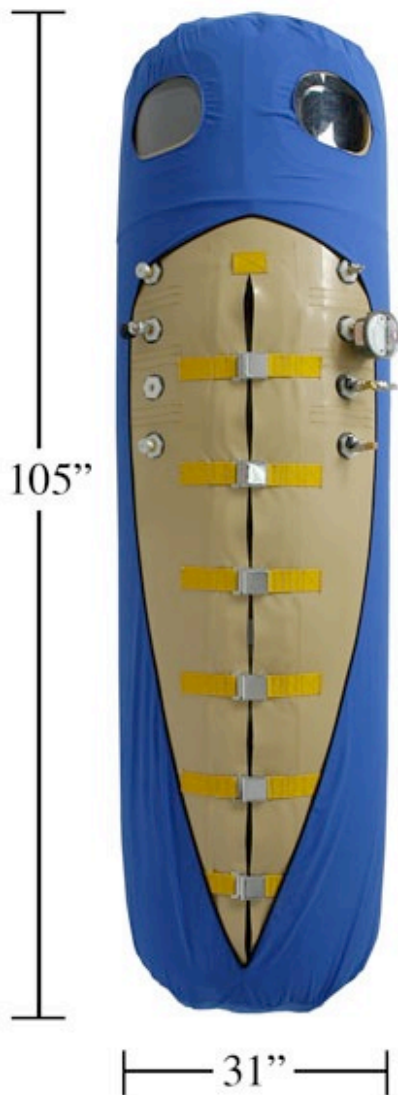

Our oversized high accuracy analog pressure gauge can be easily read from across the room. Our gauge also corrects for changes in altitude and barometric pressure insuring correct readings in any environment.



Extra long zipper and generous chamber dimensions make it possible for both a parent and child to use the chamber simultaneously. Dual blow off valves automatically vent the chamber at 4 psi preventing overpressurization.



Six stainless steel buckle closures each capable of 10,000 pounds of pull-strength are used to close the outer bag. Like the metal zipper beneath it, the buckles are easy to open and close making solo operation possible.



Double bag construction ensures superior strength and longevity. The inner bag is made from 880-weight urethane coated nylon. The outer bag is made of double urethane coated 39 oz. 2 X 2 basket weave nylon.

Our FDA 510 K approved compressor features a multiple filter system that both cools and dehumidifies the air for a more comfortable chamber environment. A custom sound baffle system decreases the sound signature of the compressor making it quiet enough to use in a medical clinic or small apartment. Compressor air quality exceeds the highest known standards regarding particulate matter size and density, carbon dioxide, carbon monoxide, ozone and 75 different VOC's.



Our 1/2 inch anodized aluminum internal frame holds the chamber in an open position even when uninflated. This makes for easier entry and exit and faster inflation time. Given the high strength of this frame, you can also use it for physical support when sliding in and out of the chamber. Quick disconnect attachments make for easy assembly. Our medical grade internal mattress is custom fitted to the chamber and is 98" long. The mattress cover unzips for easy cleaning. Closed cell anti-roll bolsters keep the chamber from rolling or rotating in place. Non-electric intercom and internal gauges also available.



To speak with a professional about Hyperbaric Oxygen Chambers Call Toll Free 1-888-824-7558